ConvaFoam™











Border

INSTRUCTIONS FOR USE / INSTRUCCIONES DE USO

Adhesive Silicone Hydrofiber™ Foam Dressing / Apósito de espuma de Hydrofiber™ con adhesivo de silicona

www.convatec.com

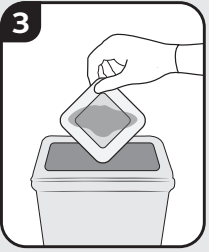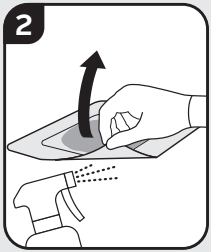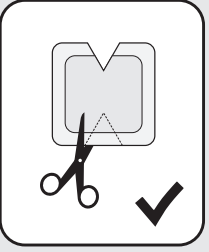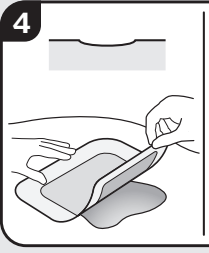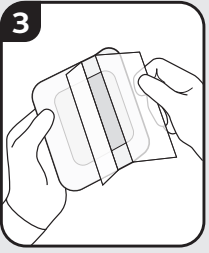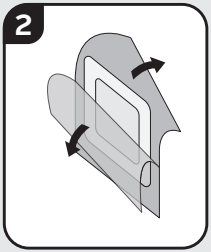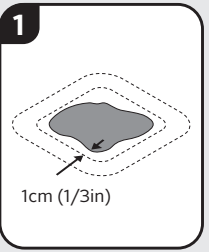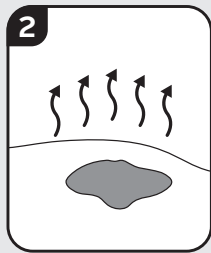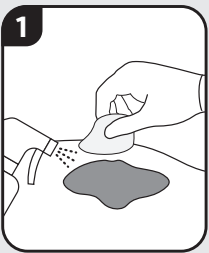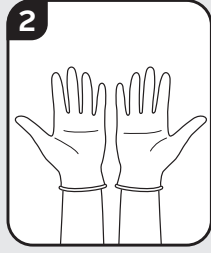PRODUCT DESCRIPTION

ConvaFoam™ Border dressings are sterile, Hydrofiber™ foam wound dressings, consisting of:

- a showerproof, breathable outer polyurethane film which provides a bacterial and viral barrier.
- a soft, multi-layered, absorbent central pad containing superabsorbent fibres, polyurethane foam and Hydrofiber™ (Sodium carboxymethylcellulose) wound contact layer. The pad absorbs and retains moisture, exudate and bacteria.
- a perforated silicone adhesive skin contact border to adhere the dressing to the skin.

INTENDED USE

ConvaFoam™ Border dressings may be used, with the consultation of a healthcare professional, for the management of exuding, non-exuding wounds and protection of intact skin.

ConvaFoam™ Border dressings may be used as a primary or secondary dressing.

TARGET PATIENT POPULATION

ConvaFoam™ Border dressings are intended to be used on patients with one of the wound types listed in the indications or those who require protection from skin breakdown.

INDICATIONS

ConvaFoam™ Border dressing is indicated for the treatment of:

- Leg ulcers
- Pressure Ulcers
- Diabetic Foot Ulcers
- Surgical Wounds
- Partial thickness burns
- Traumatic wounds

CONTRAINDICATIONS

ConvaFoam™ Border dressings should not be used on individuals who are sensitive to or who have had an allergic reaction to the dressing or components stated in the Product Description above.

ConvaFoam™ dressings are not compatible with oxidising agents such as hydrogen peroxide or hypochlorite solutions.

PRECAUTIONS AND OBSERVATIONS

ConvaFoam™ Border dressings are for single-use only. Do not reuse.

Reuse may lead to increased risk of infection, cross-contamination and delayed healing.

For highly exuding wounds, ConvaFoam™ Border may be used as a secondary dressing in conjunction with a primary dressing, such as AQUACEL™ silver and non-silver Hydrofiber™ products.

ConvaFoam™ Border dressings should not be submerged during bathing or swimming.

Consult a healthcare professional if any of the following are observed:

- Irritation (reddening, inflammation)
- Maceration (whitening of skin)
- Hyper-granulation (excess tissue formation)
- Sensitivity or allergic reaction
- Signs of infection (increased pain, bleeding, warmth/redness of surrounding tissue, wound exudate)
- Changes in wound colour and/or odour
- Increased wound size
- Wound not showing signs of healing

After use, this product may be a potential biohazard.

Handle and dispose in accordance with accepted medical practice and applicable local laws and regulations.

If during the use of this device or as a result of its use a serious incident has occurred, please report it to the manufacturer at www.convatec.com and to the competent authority of the Member State in which the user and/or patient is established.

CONTACT AND DURATION OF USE

ConvaFoam™ Border dressings can be worn for up to 7 days. Dressings should be changed earlier if clinically indicated. It is safe to continue using ConvaFoam™ Border dressings for as long as clinically indicated.

DIRECTIONS FOR USE

1. Storage

- Store at room temperature.
- Protect from light.
- Keep dry.
- Ensure the dressing pouch is not damaged or opened prior to use to guarantee sterility. If the dressing or pouch is damaged or opened, do not use and dispose according to local regulations.
- Do not use after expiry date.

2. Wound Site Preparation and Skin Cleansing:

- Ensure the skin surrounding the wound is free from creams and moisturisers which could prevent the dressing from sticking.
- Before applying the dressing cleanse with a suitable wound cleanser and dry the surrounding skin.
- If a moisture barrier is required, select a spray or liquid that will dry before dressing application.

3. Dressing Preparation and Application:

- When applying to a wound choose an appropriate dressing size where the central absorbent pad is at least 1cm (1/3 inch) larger than the wound.
- When applying to intact skin choose a dressing which covers the entire vulnerable area.
- Remove the dressing from the sterile pouch.
- Remove the release liner. Try to avoid touching the adhesive surface.
- Hold the dressing over the wound or area of skin to be

- protected and line up the centre of the dressing with the centre of the wound or area of skin to be protected.
- Ensure the dressing is not stretched upon application. Smooth down the adhesive border. Ensure joints are sufficiently flexed to allow patient mobility.

If required, dressings may be cut to fit. Use additional tape or film dressing to create a seal.

4. Dressing Removal:

- The dressing should be changed when clinically indicated, when saturated with wound fluid, if it leaks, or if it becomes damaged or soiled.
- If removing the dressing is difficult, the dressing should be soaked with sterile saline or water and removed slowly.
- To remove the dressing, press down gently on the skin immediately surrounding the dressing and carefully lift one corner of the dressing. Continue lifting the dressing until all edges are free. Carefully lift away the dressing and discard as per local clinical protocols.
- During protective care only, the dressing corner may be lifted for inspection of skin and resealed. The dressing must be changed if skin is broken.

If further information or guidance is needed, or to report an adverse event, please contact ConvaTec, contact details can be found at www.convatec.com or use the contact details at the end of this document.

Made in UK

© 2022 ConvaTec Inc.

®/™ indicates a trademark of ConvaTec Inc.

Caution: Federal law restricts this device to sale by or on the order of a physician

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®/™ indica una marca registrada de ConvaTec Inc.
ConvaTec Limited